ELECTRIC TOOTHBRUSH

THANK YOU FOR CHOOSING US

Smart reminder | Cultivate the habit of scientific tooth brushing Whole body washing | Safe to use, easy to clean Sonic vibration | Improve the efficiency of brushing, reduce bleeding of your gums DuPont bristles | International brand, safe using Long-term use | Clean dead corners between teeth, brighten and clean teeth Keep using | Avoid oral diseases and bring you a confident smile

We are always committed to let every consumer can easily enjoy high-quality

oral care in a simple way. Experience! Effectively clean stubborn tooth stains such as dental plaque and tobacco stains, tea stains, coffee stains, etc.

Prevent gingivitis, relieve red and swollen gums, and improve the oral hygiene of orthodontic patients. For your convenience, Please read the instructions for use of this product in detail.

Thank you very much for choosing this product. We have been focusing on the production and development of electric toothbrush related products. Not only to make good and designable products. What we really want to do is to guide the electric toothbrush correctly consumers' attitudes towards teeth and oral care. Through our efforts, let you and your family stay far

IMPORTANT SAFETY Notice

1. Do not immerse or place the charging cable in water or any liquid. If you find that the charging cable has fallen into the water ,please unplug the plug immediately. Do not reach out to remove the charging cable that has fallen into the water. After cleaning, make sure the charging cable is finished After it is completely dry, turn on the power again.

 This product contains a built-in rechargeable battery, please do not throw it into a fire or charge, use or put it in a high temperature place. 3. Please keep away from strong electromagnetic field when charging this product.

Warning

1. Before connecting the power supply, please check whether the local power supply voltage is consistent with the voltage marked on the product.
 This product contains parts that cannot be repaired by the user (such as:

charging cable, toothbrush handle, wire, etc.)
If any product is damaged, please replace the access
corresponding model in time to avoid danger.

3. Please place it out of the reach of children and babies, and cannot be used as a toy.

4. If you have had oral or gum related surgery within two months, please

consult your dentist before using this product 5. This product is only suitable for cleaning teeth, gums and tongue, please

This product is only suitable on cleaning teeth, goins and congue, product of once use it for other purposes.
 If the toothpaste contains hydrogen peroxide, baking soda or bicarbonate, please wash your brush head thoroughly with water after use it
 Do not use this product in a bath or bath. If you do not use it for a long

time, please unplug the power plug to avoid fire or danger.

CHARGE DESCRIPTION

1. Please check the power supply voltage and power socket $% \left(x\right) =\left(x\right) +\left(x\right) +$ voltage adaptation

2. Before plug the appliance into the power supply under safe conditions, please make sure that the charging cable is dry and the corners of the

Please read the product safety tips before using, and keep it carefully. 01

ELECTRIC TOOTHBRUSH

charging plug are free from water stains and damage.

3. Please power on the electric toothbrush before first using it. In case of low power or no power of the toothbrush, please charge it., After start-up, put the toothbrush into the connected power supply to charge, and the toothbrush stops vibrating and enters charging mode.

Battery power capacity	Power indicator status	LCD Display
Charging	LCD screen displays the current battery level	Changing 45%
Charge complete	LCD screen displays the current battery level	Омизия 100%
Power less than 10%	Low power prompt, red lightning slow flashing animation, start-up animation appears after 3 seconds, normal function settings can be performed by operating the function key within 8 seconds, and it will automatically shut down if there is no operation within 8 seconds	Low power Start-up prompt animation
Power less than 3%	Power on prompt animation, red lightning quick flashing warning power run out. It cannot enter the functional operation interface in the time. Automatic shutdown lock after 8 seconds. Only power more than 3%, it can be normally turn on and use	Low Promote

When the battery is in low power protection mode, the red lightning icon on the LCD screen flashes quickly, and the toothbrush can only opens but motor doesn't work.

1. We elaborately designed the battery indicator for you. Please pay attention to the status of the battery indicator during use. Please charge it in time when the battery is low.

If you have not used the electric toothbrush for a long time, please charge the toothbrush first when you use it again. And to ensure safe and correct charging, which is conducive to improving the service life of the product.

Please do not dispose of the product randomly after the end of its service life. The built-in battery may contain ingredients harmful to the environment. Please dispose of the product in the battery disposal area of your area and ensure that the internal power of the battery is exhausted.

Please read the product safety tips before using, and keep it carefully. 02

I ELECTRIC TOOTHBRUSH

OPERATION AND FUNCTION DESCRIPTION

1. Mode switching: touch or key to turn on and enter the boot interface. Short press the mode function selection key to enter the mode selection. (Short press the switch button to enter this mode to work). Sequential function mode: During the toothbrush working process, short
press the mode function selection key to switch the function modes in cycles
(Cool Clean---Sensitive and Soothing--- Deep Cleaning---Whitening and Polishing---Comfortable Massage---DIY ---setting)



O.

DIY

3.On&Off

Comfort massage

a. Touch the trigger area on both sides , the electric toothbrush will light up $% \left\{ 1,2,\ldots,n\right\}$ and enter the startup screen, short press the switch button within 8 seconds to enter the working mode test state when it was turned off last time, and automatically enter the shutdown state if there is no operation within 8

(D)

b. Short press the power button or the function mode button, the electric but on the less till epoten but on the factor of the factor to be but on, the factor to though the but on and enter the start up screen, short press the mode button to work in the last shutdown mode, short press the switch button, the electric toothbrush enters this mode, press the switch button repeatedly, the electric toothbrush appears 8 seconds after the shutdown



(·:) 0:00

Please read the product safety tips before using, and keep it carefully. 03

ELECTRIC TOOTHBRUSH

4. Setting function: After entering the setting interface, press the confirm key, and then press the mode selection key to enter the system function setting. At this time, you can perform the system language, restore the

factory and other function settings. Language setting: This product supports bilingual switching between Chinese and English. Setting interface-language selection for setting $\label{lem:Restore} \textit{Restore factory settings: Settings interface-restore factory settings interface}$



Packing List

Charger Cable





PRODUCT SPECIFICATIONS & INFORMATION

Model	SN-P1	
Color	White, black (specific to the actual product)	
Body		
Rated Voltage/Power	3.7V 2W	
Duration of service	Approximately 2 hours (calculated in deep cleaning mode)	
Standard	GB4706.1-2005, GB4706.59-2008	
Frequency	36000RPM/min (calculated by deep clearing mode)	
Battery	Lithium battery	
Charging time	about 3 hours	
Product size	27.5*245.5mm	

Sonic electric toothbrush

Note: If this product is changed, please pay attention to the official information without prior notice.

to change the brushing area. 6. Automatic memory function, it will be the previous shutdown memory gear mode operation after starting up. Toothbrush head -LCD display ·

5. Intelligent brushing reminder, automatically shut down after 2 minutes; In the working status, there is 0.5 second short stop every 30 seconds, remind

ELECTRIC TOOTHBRUSH



regional reminde

USB charging







USB charging (TYPE-C interface)



IPX7 waterproof






DIY selection

Please read the product safety tips before using, and keep it carefully. 05

| ELECTRIC TOOTHBRUSH


BRUSHING TIME

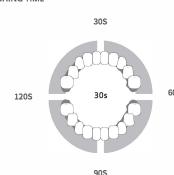

- $\bullet\,$ The dentist recommends that you should brush your teeth for 2 minutes
- brush their teeth, which has achieved a scientific cleaning effect.
- . After each use, the toothbrush will remember the selected mode, when the next time turn on, the memory pattern begins to execute.

QUICK USING GUIDE >

Step 1: Before using the toothbrush for the first time, please start the toothbrush. If there is any remaining battery, please use it up and recharge it. You can charge it as you use it later, but it is also recommended that you fully discharge the battery once a month to extend the service life of the lithium

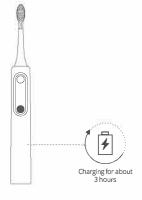


Please read the product safety tips before using, and keep it carefully. $\hfill \hspace{0.5cm}$ 06

| ELECTRIC TOOTHBRUSH


Step 2: Install the brush head firmly on the brush head shaft. There is a gap of about 1mm between the brush head and the handle to avoid resonance tha affects the experience during the brushing



Step 3: Squeeze an appropriate amount of toothpaste on the toothbrush

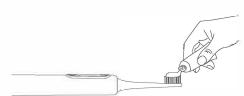


Step 4: Place the toothbrush gently at the junction of the teeth and gums

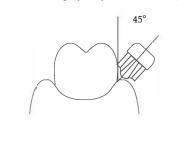


Please read the product safety tips before using, and keep it carefully. $\hfill 07$

I ELECTRIC TOOTHBRUSH


Step 5: Press the power button to start cleaning your teeth and clean all parts.

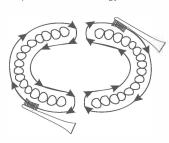

Step 6: To prevent the toothpaste foam from splashing, please turn off the toothbrush before taking it out of your mouth

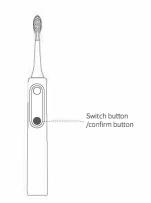

Step 7: After brushing the teeth, rinse the toothbrush head with running water, shake off the water on the brush head, and wipe the handle clean $\frac{1}{2}$

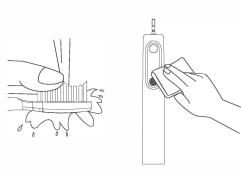

Please read the product safety tips before using, and keep it carefully. 08

| ELECTRIC TOOTHBRUSH

Step 8: Put the electric toothbrush upright on a dry table, remember to put on the dust cover of the brush head.

Please read the product safety tips before using, and keep it carefully. 04



Storage

1. If you are not about to use the electric toothbrush for a longtime, please unplug the charger from the power socket, clean it and store it in a cool place., Dry and away from direct sunlight.

2. The electric toothbrush is suitable for indoor use, and is suitable for storage in a ventilated and dry place at 0-35 degrees Celsius. 3. Regularly clean the power plug of the charging cable to prevent the $\,$ accumulation of dust and dirt and reducing the insulation effect which can lead to overheating or electric shock causing unnecessary loss. 4. Do not place metal objects, such as coins or paper clips, near the toothbrush.

5. If you don't use it for a long time (recommended not to charge for more than 12 hours), please unplug the power supply to avoid fire or danger. 6. This product contains lithium batteries. Please store them away from fire and heat sources. Do not charge, transport, use and place in high emperature environments. 7. If the toothpaste contains hydrogen peroxide, baking soda or carbonate

(whitening toothpaste), please clean it thoroughly with water and place it in a cool and dry place after brushing your teeth

Sonic Electric Toothbrush

Instructions

Please read the product safety tips before using, and keep it carefully. $$ 09